



University Of AL-Qadisiyah

DEPARTMENT OF MEDICAL PHYSICS

PhyDCM Application Helps

Left panel

Patient Name:

Enter patient name

Patient Age:

Yo

Gender:

Male

Scan Type:

MRI

Model Status:


PhyDCM (Default)

Import Medical File

Import DICOM Folder

Choose

Start Diagnosis

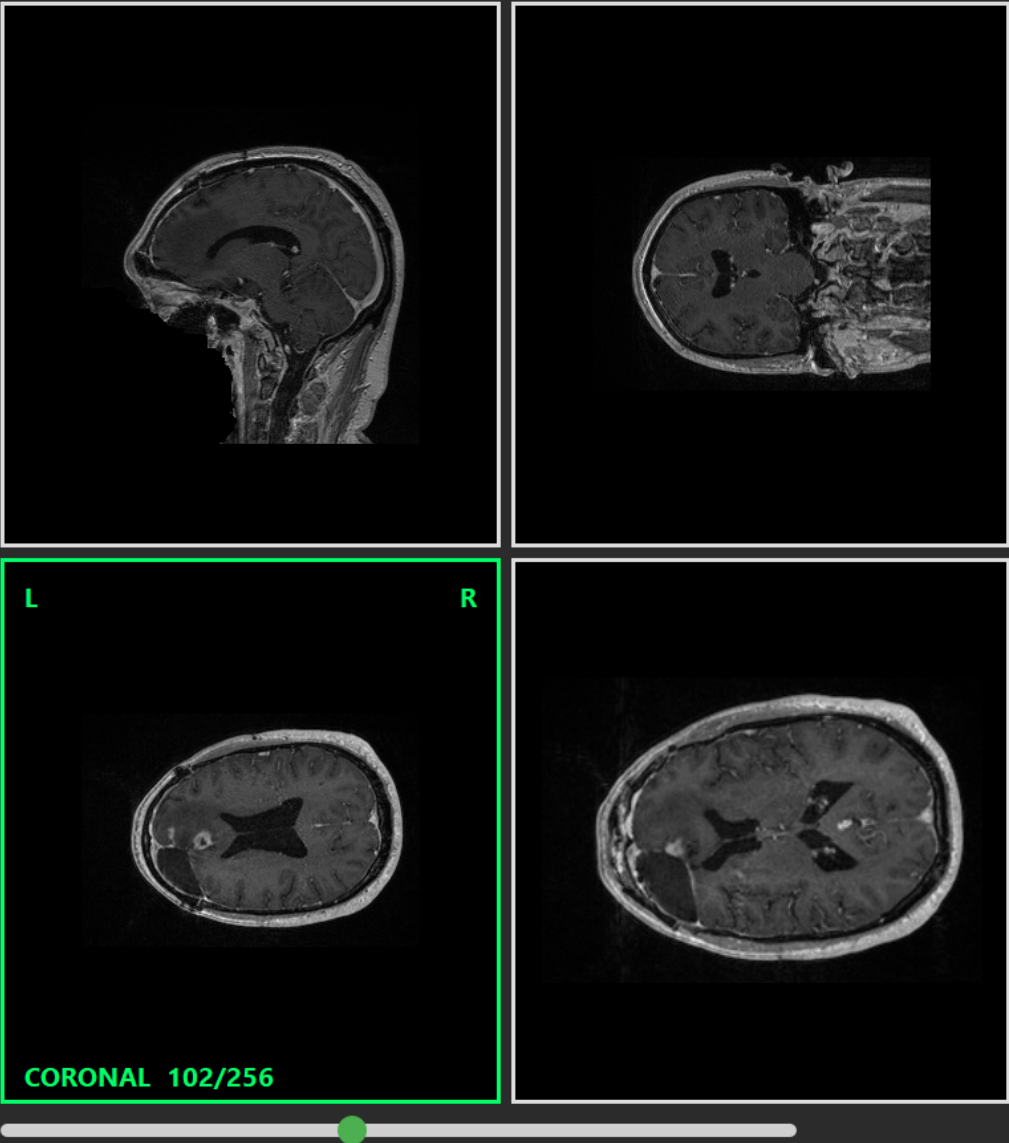

	Sequence Count	ID	Patient Name	Age	Gender
1	15	215021M0111	Mohammed Hadi	21	Male
2	176	112625G01G1	ggggg	25	Male

Export CSV/JSON

Add Patient

Delete Selected

Left panel

It includes a patient name field:

The name must be between three and fifteen characters long and must be in English only, without any symbols, numbers, or other languages.

The age field will later appear in the table in the right-hand box.

The gender selection field is limited to two options: male/female.

The scan type selection field, currently in this version of PhyDCM, includes three devices: CT/MRI/PET. Below this field, the PhyDCM library, trained on these models, is displayed in orange.

The upload (folder) button opens 4D animated images on the four screens within the middle box, while the upload (file) button opens a single image on the fourth screen.

The "Choose" button selects the image from the folder suitable for diagnosis (it only works when an image folder is uploaded), and the "Start diagnosis" button begins the diagnostic process.

File Edit View Import Info

Patient Name:

Enter patient name

Patient Age:

Yo

Gender:

Male

Scan Type:

MRI

Model Status:

PhyDCM (Default)



Import Medical File

Import DICOM Folder

Choose



Start Diagnosis

Not Available

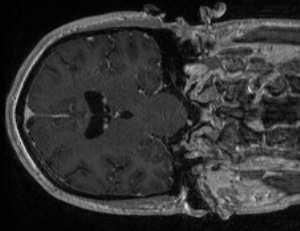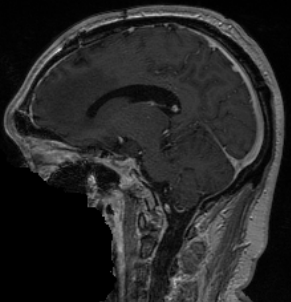

L

R

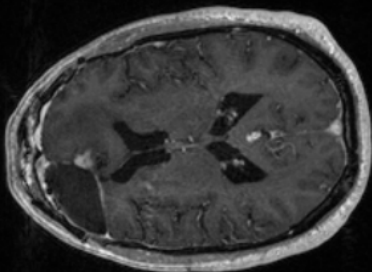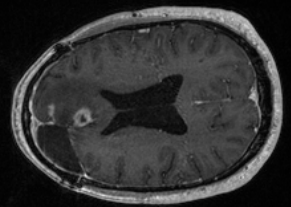

CORONAL 102/256

Sequence Count

ID

Patient Name

Age

Gender

1

15

215021M0111

Mohammed Hadi

21

Male

2

176

112625G01G1

ggggg

25

Male

Export CSV/JSON

Add Patient

Delete Selected

Dark mode is enabled

File	Edit	View	Import	Info
New Add	Ctrl+N			
New Window	Ctrl+Shift+N			
Open File	Ctrl+O			
Save	Ctrl+S			
Save As...	Ctrl+Shift+S			
Share				
Auto Save				
Clear File	Ctrl+K, F			
Close Window	Ctrl+F4			

File	Edit	View	Import	Info
	Undo (Not Available)	Ctrl+Z		
	Redo (Not Available)	Ctrl+Y		
	Copy	Ctrl+C		
	Paste	Ctrl+V		
	Zoom In	Ctrl++		
	Zoom Out	Ctrl+-		

Info
About Us
Help

Note that the application only accepts uploads in medical formats such as DCM or NII, and it also counts the number of images within the uploaded folder.

Other features not currently available in this version include the import function, which uploads a library or training files that can be used for diagnosis, while the primary library within the application is PhyDCM.

Address Bar

File	Edit	View	Import	Info
		Processing (Not Available)		
Patient		Terminal	Ctrl+T	
Enter		History		
Patient		Segmentation (Not Available)		
Go		Light/Dark Pattern	Ctrl+D	

Patient Name:

Enter patient name

Patient Age:

Yo

Gender:

Male

Scan Type:

MRI

Model Status:


PhyDCM (Default)

Import Medical File

Import DICOM Folder

Choose

Start Diagnosis

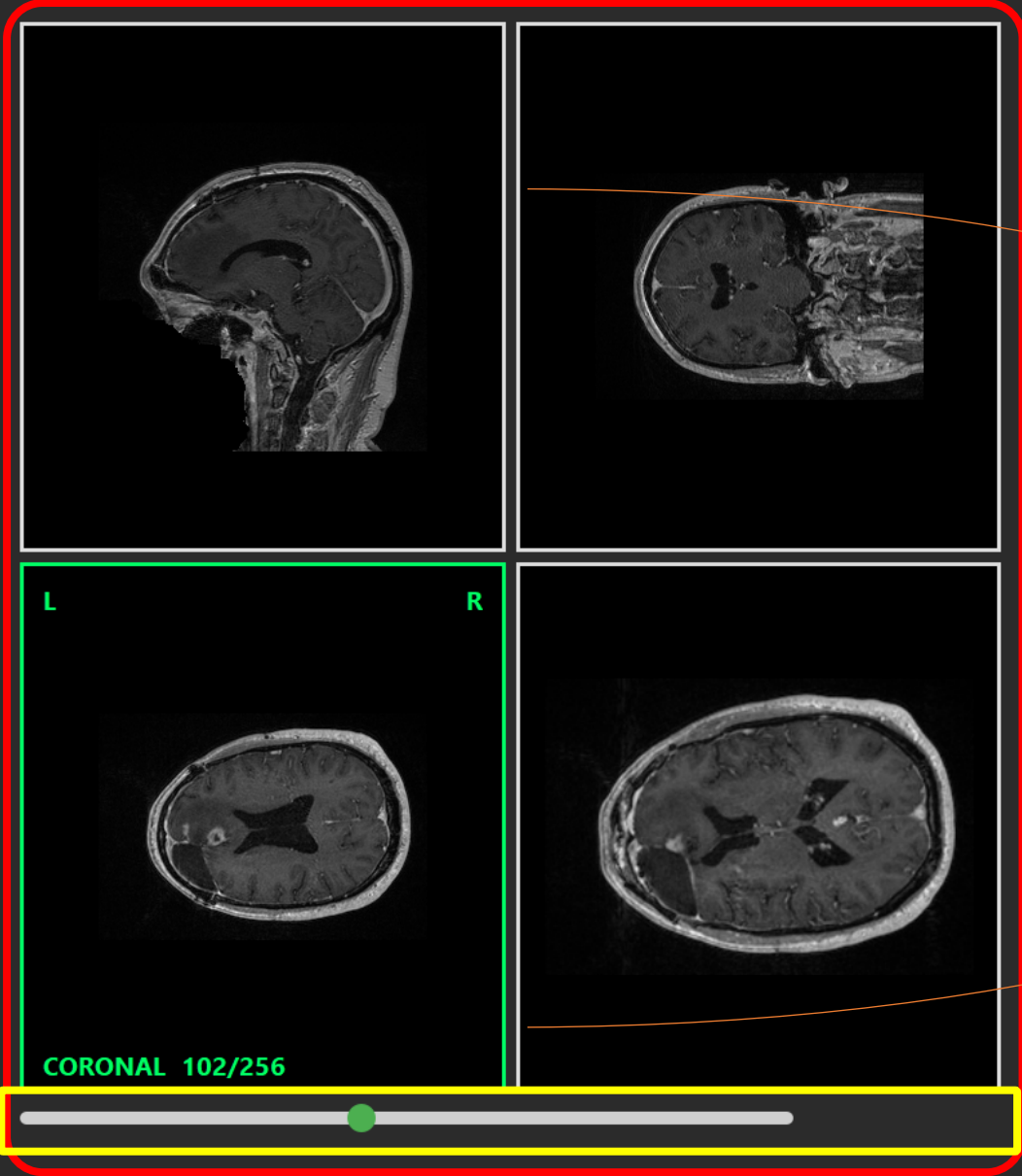

	Sequence Count	ID	Patient Name	Age	Gender
1	15	215021M0111	Mohammed Hadi	21	Male
2	176	112625G01G1	ggggg	25	Male

4 Screens,
AXIAL/SAGITTAL/CORONAL

The main screen
scroll bar has a
green border

Right panel

Patient Name:

Enter patient name

Patient Age:

Yo

Gender:

Male

Scan Type:

MRI

Model Status:


PhyDCM (Default)

Import Medical File

Import DICOM Folder

Choose

Start Diagnosis

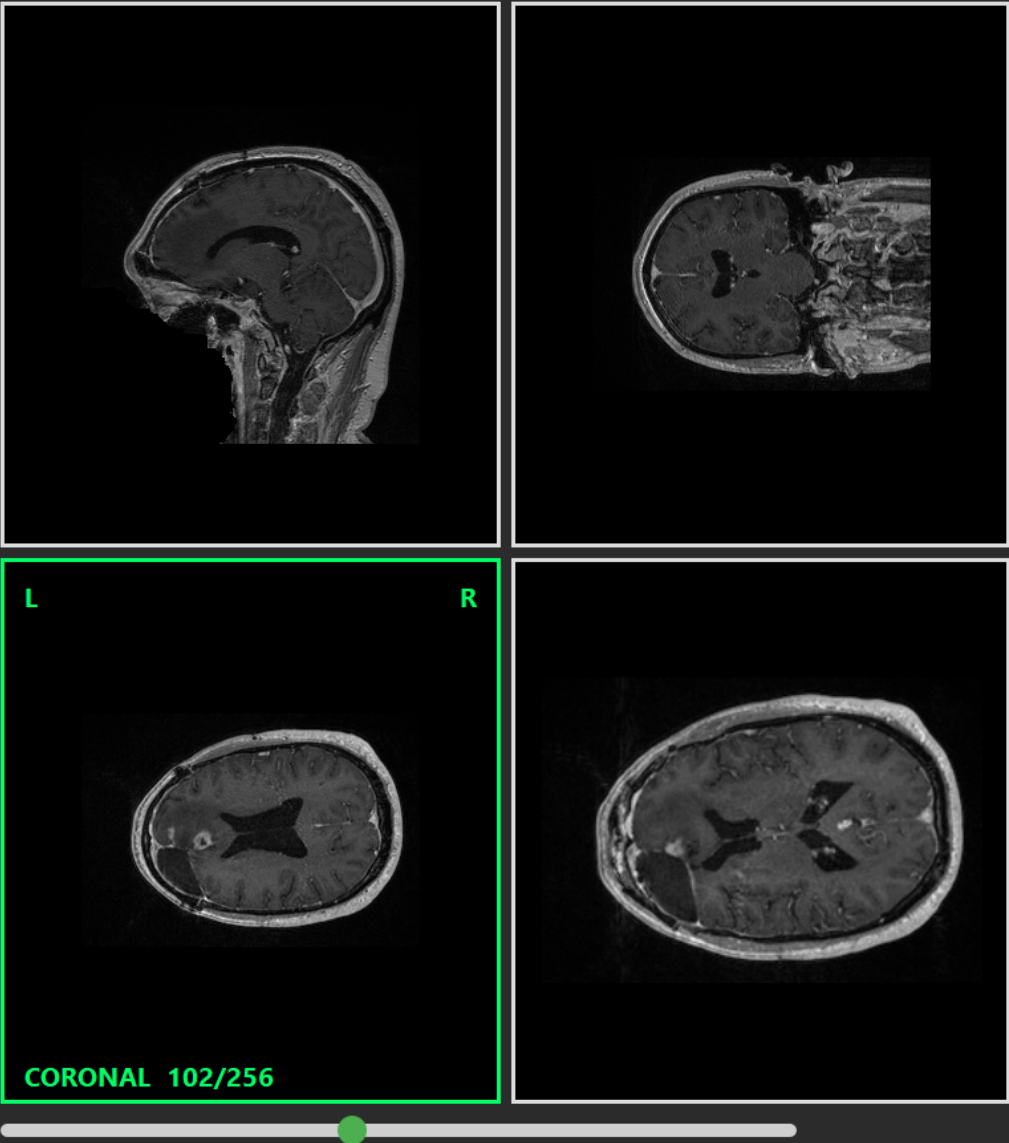

	Sequence Count	ID	Patient Name	Age	Gender
1	15	215021M0111	Mohammed Hadi	21	Male
2	176	112625G01G1	ggggg	25	Male

◀

▶

Export CSV/JSON

Add Patient

Delete Selected

Right panel

After the diagnostic process begins, all the information and image features will be calculated, in terms of entering the fields of names, the patient's ID, gender, age, type of examination, image size and type, and even its dimensions.

Then it comes to other details such as brightness and sharpness, and then the AI-powered diagnostic process appears based on the virtual library PhyDCM.

The number of seconds of diagnosis appears in its last fields, and then it can be saved as CSV, Excel, or Json (it appears in the history box).

Export CSV/JSON

Add Patient

Delete Selected

HISTORY Box

Search ID: Type Patient ID...

Newest First ▾

Refresh

112625G01G1

215021M01I1

Newest First

Oldest First

Name A-Z

Name Z-A

Load Selected

Open HISTORY Folder

Patient Name:

Enter patient name

Patient Age:

To

Gender:

Male

Scan Type:

MRI

Model Status:


PhyDCM (Default)

Import Medical File

Import DICOM Folder

Choose

Start Diagnosis

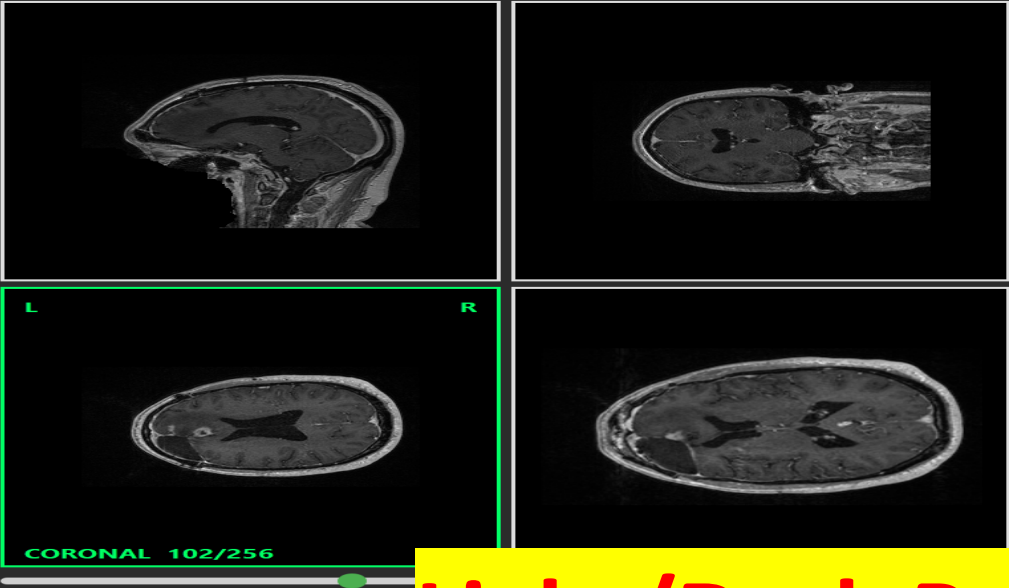

	Sequence Count	ID	Patient Name	Age	Gender
1	15	215021M01I1	Mohammed Hadi	21	Male
2	176	112625G01G1	ggggg	25	Male

Export CSV/JSON

Add Patient

Delete Selected

Light/Dark Pattern

Patient Name:

Enter patient name

Patient Age:

18

Gender:

Male

Scan Type:

MRI

Model Status:


PhyDCM (Default)

Import Medical File

Import DICOM Folder

Choose

Start Diagnosis

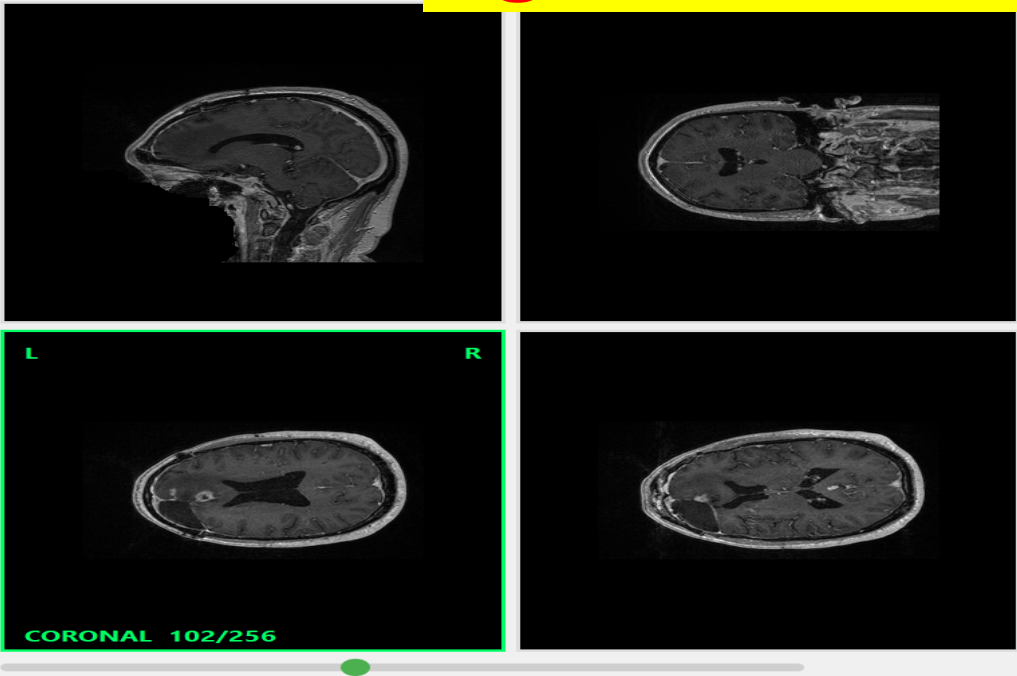

	Sequence Count	ID	Patient Name	Age	Gender
1	15	215021M01I1	Mohammed Hadi	21	Male
2	176	112625G01G1	ggggg	25	Male

Export CSV/JSON

Add Patient

Delete Selected